HF- Bipolare Elektroden

HF bipolar electrodes

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Seite 3

Operating Manual

Page 15

DEUTSCH


ENGLISH

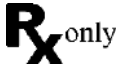Inhalt

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Web: www.HEBUmedical.de

Contents

1	Symbol descriptions	16
2	Introduction	17
3	Scope	17
3.1	Intended Use	17
3.2	Contraindication	17
3.3	Maximum electrical capacity	17
4	Instructions for Use	17
4.1	Connection to HF unit in bipolar mode	17
5	Warnings	18
6	Safety Remarks	19
7	Preparation	20
7.1	Reutilization restrictions	20
7.2	Information on instrument preparation	20
7.3	Preparation at the place of use	20
7.4	Ultrasound bath (optional)	20
7.5	Manual cleaning	21
7.6	Mechanical cleaning	21
7.7	Drying	23
8	Maintenance, inspection	23
8.1	Functional check	23
9	Sterilization	24
9.1	Packaging.....	24
10	Lifetime	24
11	Storage	25
12	Warranty/Repair	25
13	Servicing and manufacturer address	25

1 Symbol descriptions

Symbol	Definition
	CE-labelling
	Attention
	Validated Parameters
	Manufacturer
	Lot-description
	Reference code
	Medical device / FDA Prescription device
	Medical device
	Non sterile
	Keep away from sunlight
	Dry storage required
	(Electronic) instruction for use
	No intervention on patients with cardiac pacemakers.
	Temperature

2 Introduction

Our products are exclusively intended for professional use by appropriately trained and qualified personnel and may only be acquired by them.

By purchasing this instrument, you are now the owner of a high-quality product whose use and correct handling are described in the following. In order to minimize possible risks to patients and users, please observe these instructions carefully. Use, disinfection, cleaning and sterilization may only be performed by suitably trained specialist personnel.

3 Scope

The use of HF bipolar electrodes is used in HF surgery for coagulation, volume reduction and soft tissue cutting.

3.1 Intended Use

A bipolar or monopolar electrical conductor (blade or tip) connected to a mating electrical handpiece for establishing an active electrical connection between the output terminals of an electrosurgical generator and a patient. The delivered energy is used to perform lesions (cutting and / or cauterizing) at a desired location on the patient. It is commonly used in various medical disciplines. It is a reusable product.

3.2 Contraindication

The HF bipolar electrodes should not be used in conjunction with the central nervous or circulatory system. Tube ligation is not a safe application.

3.3 Maximum electrical capacity

Item no.	U_{\max}
HF9596-01 to HF9697-03	0,5 kVp
HF9698-01 to HF9699-03	1,0 kVp

4 Instructions for Use

Thoroughly examine all products for visible irregularities. Do not use cables with a brittle or defective insulation. In case of doubt check operation by means of a suitable testing device.

4.1 Connection to HF unit in bipolar mode

In bipolar HF surgery the coagulation is restricted to the area between the two electrodes of the bipolar instrument. Consequently, a neutral electrode is necessary. Otherwise, follow the recommendations regarding patient preparation in monopolar operating mode and the corresponding instructions for the HF unit. When using the bipolar “Auto start” mode, the cable used should not be longer than 4.5 m. Alternatively, use a twin cable with max. length of 5m.

5 Warnings

The products may only be used in combination with HEBU accessories and only by clinically instructed and qualified personnel. Electrodes must not be twisted, as this may cause damage to the material and the electrode might break. Power output of the HF unit: Always select the lowest possible output of the HF unit for any procedure. If the coagulation capacity of the electrode is below normal, do not increase the output of the HF unit without a thorough preliminary check. In particular, the following must be checked: correct attachment of all HF cables and plugs, correct activation of current by the finger - or footswitches, undamaged insulation of the HF cables and instrument, and check that the distal end of the electrode is dean and undamaged. Before operating, check that the signal transmission of the monitors works without interference (e.g. without noise); this check is made by pressing the yellow button (CUT) and the blue button (COAG).

Always unplug the HF-cable by holding on the connector. **Never pull the cable itself**, as this may lead to defects. Cables with defective/brittle insulation or with a broken electrical wire (due to severe bending) might cause a fire. Active HF instruments and HF handles must not be placed on the patient in order to avoid injuries in case of involuntary HF activation or activated "Auto Start" mode.

For cleaning the active electrode, you must deactivate the "Auto Start" mode (only in bipolar applications) or disconnect the instrument from the generator. Directly after activation the tips could be hot. This may lead to burns.

Endogenous burns are those caused by high current density in the patient's tissue. Possible reasons are: The patient is inadvertently positioned in contact with electrically conductive parts. Direct contact between skin areas and HF cables may lead to electrical capacitance which in turn may cause burns.

Exogenous burns are those caused by heat of ignited fluids or gases. They may also be caused by explosions. Possible reasons include: Ignition of skin cleaning agents and disinfectants, ignition or narcotic gases etc.

Cardiac pacemakers may be damaged by electrosurgical current. Consult a cardiologist before operating. Never use electrosurgical current on outpatients fitted with pacemakers.

6 Safety Remarks

	Check the isolation on the electrode and the shaft extensions for damage before each use.
	Do not rebend needles, loops or knives. Spatulas can be bent.
	If damaged, do not use these products!
	Do not use in the presence of combustible or explosive materials!
	Never lay instruments on the patient or in his/her direct vicinity!
	HEBUmedical expressly warns against modifying the product. Any modification exempts HEBUmedical in from any and all liability. .
	Protect these products from any form of mechanical damage! Do not throw! Do not: use force!
	Defective HF cables or electrode handles may not be repaired. Always replace a defective part with a new one.
	A visual inspection of the instrument for damage and contamination must be carried out before each use!
	Avoid mechanical overstressing of the instrument beyond the structural design, this can lead to breakage and deformation!
	If the products are used on patients with transmissible spongiform encephalopathy or HIV infection, we decline any responsibility for their reuse.
	Observe manufacturer's instructions for electrical load capacity!
	Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

7 Preparation

The person in charge of preparatory treatment is responsible for ensuring that the treatment is duly carried out using the relevant equipment, materials and personnel in the treatment facility and so achieves the desired result. This necessitates validation and routine monitoring of the process used. We urge you to take note of the national regulations dealing with instrument preparation.

7.1 Reutilization restrictions

Frequent repeat preparatory treatment has minimal effects on the product. The end of the product life is normally determined by wear and damage due to use

7.2 Information on instrument preparation

- Use cleaning and/or disinfection agents with a pH-value within 9-10. Please observe manufacturer instructions regarding dosage, exposure time and renewal of solutions.
- Do **not** use hard brushes (e.g. metal brushes or metal sponges) or coarse abrasive cleaners.
- Never leave instruments in cleaning or disinfection agents for longer than the specified time.
- Only used demineralized water for rinsing.
- Rinse and dry carefully through channels and pipes.
- Sensitive instruments must be cleaned in a storage or clamping fixture.
- Observe manufacturer instructions of cleaning – and sterilizing equipment.

7.3 Preparation at the place of use

Directly after using remove coarse dirt of products and rinse out the working cannulas. Do not use fixing agents or hot water (> 40°C), as this results in residues becoming fixed and can affect the success of the subsequent cleaning operation

Dismantle and/or open products as far as possible. Within short time after use clean the products for reducing a drying of the residues.

This enables an easier cleaning. If products come into contact with corroding medicines or cleaning agents, wash these up with water immediately after use.

Longer drying times, e.g. for dry disposal are not validated and therefore not recommended.

 The drying time during validation was 1 hour.

7.4 Ultrasound bath (optional)

All products must be opened, dismantled and any cavities rinsed through.

Place products in the screen basket in such a way that overlaps and contact between products are avoided. Add cleaning agent to the water and adjust the temperature of the solution in line with the cleaning agent manufacturer's instructions.

The cleaning in the ultrasound bath should be at 35-40 kHz, 5 minutes at least.

 To validate cleaning in an ultrasonic bath, the test items were ultrasonically treated in Neodisher Mediclean forte 0,5 % for 5 minutes.

Subsequently rinse products including all cavities before cleaning and disinfection.

Medicine products which are possessing a bad ultrasound transmission, e.g. soft material are not usable for the ultrasound bath.

7.5 Manual cleaning

Since mechanical processes can be standardized, reproduced and therefore validated, mechanical cleaning/disinfection should be preferred to manual processes.

Manual cleaning and disinfection process is not validated and therefore needs to be validated additionally by the end user.

7.6 Mechanical cleaning

On the basis of international standards (EN ISO 15883) and national directives, only validated machine cleaning and disinfection methods may be used. For the mechanical cleaning we recommend a standard program.


Only completely demineralized water should be used for cleaning, neutralisation and rinsing, in accordance with the „Guidance Compiled by the DGKH (Germany Society for Hospital Hygiene), DGSV (German Society for Sterile Supply) and AKI (Working Group on Instrument Reprocessing) for the Validation and Routine Monitoring of Automated Cleaning and Thermal Disinfection Processes for Medical Devices as well as Advice on Selecting Washer-Disinfectors“ (which refers to DIN EN ISO 15883-1 Point 6.4.2)

Flexible (complex) products with invisible surfaces have to be pre-cleaned manually before mechanical cleaning.



We recommend for all push shafts, -and pipe shafts instruments and instruments whose surfaces are on top of each others during the cleaning (f.e. bone forceps and gouge forceps) a manually pre-cleaning for an optimal cleaning result without residues

Observe the following by loading:

- Place the dismantled/opened instruments securely in the tray.
- Products with openings and gaps have to be faced down with the opened side so that they can be cleaned and no water of the cleaning process is collecting inside them. If available use balanced devices for rinsing

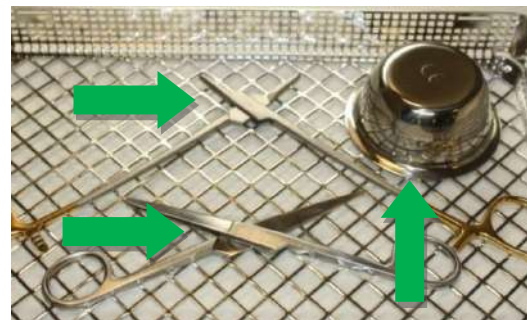- Do not overload trays, avoid creating any overlaps.

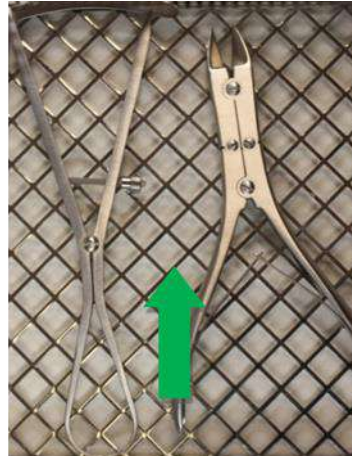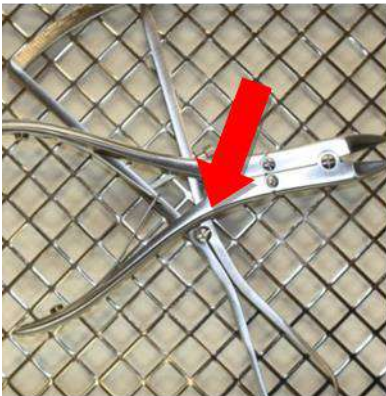Preliminary rinsing (cold, if applicable fully demineralized water without additives) is followed by chemical. The chemical cleaning should take place at **40°C -60°C** for at **least 5 minutes**.

We recommend products with a **pH-value within 9-10**, e.g. Neodisher MediClean forte from Dr. Weigert. The cleaning agents used should be selected depending on the material and properties of the instruments and in accordance to national regulations: If there is a high chloride concentration in the water, pitting and tension crack corrosion can occur on the products. The occurrence of this type of corrosion is minimized by using alkaline cleaning agents and demineralized water. By adding an acid-based neutralization medium, the rinsing off of alkaline cleaning agent residues is facilitated during the first intermediate rinsing process (warm or cold water).

In order to prevent the formation of deposits, it is advisable to use neutral cleaners where the water quality is unfavourable. After the second intermediate rinsing process, thermal disinfection takes place.

The thermal disinfection should take place at temperatures of between **80 and 95°C**, with an **exposure time as outlined in EN ISO 15883**.

After the finishing of the program take the good out of the machine because corrosion can arise if the instrument remains in the machine.

 Parameters used for the validation of preparation	
Pre-rinsing	1 minute with cold tap water
Cleaning	Temperature: 55 °C
	Soaking Time: 5 minutes (worst case)
	Neodischer Mediclean forte 0,4% (worst case)
Neutralization	Temperature: cold DI water
	Soaking Time: 2 minutes
	Neodisher Z 0,1%
Post-rinsing	2 minutes with cold DI water
Disinfection	Temperature: 90 °C (A ₀ 3000)
	Soaking Time: 5 minutes

7.7 Drying

Ensure adequate drying by the cleaning and disinfection device or using other suitable measures.

V Drying was omitted in the validation (worst case condition).

8 Maintenance, inspection

After cooling to room temperature, the products must be visually inspected for protein residues and other contamination, paying particular attention to cavities, blocks, inclusions, pipes, and other inaccessible areas. Instruments which are not free of residues must be returned for a complete retreatment process.

To ensure that the products can be used according to their intended purpose after reprocessing, it is necessary to carry out a functional test after cleaning, disinfection and drying of the visual inspection and care measures. Carry out functional checks described in point 8.1.

Products with stains, which are blunt, bent, no longer function or which are otherwise damaged must be segregated!

To help identify faulty products that need to be sorted out, we recommend the brochure "Instrument Reprocessing" from the Working Group "Instrumenten Aufbereitung". This includes in particular Chapter 8 "Checks and Care" and Chapter 12 "Surface Changes: Deposits, Discoloration, Corrosion, Aging, Swelling and Stress Cracks".

8.1 Functional check

A newly purchased product must be subjected to a thorough visual and function check after its delivery and before each use. If a new product has externally visible defects or does not operate according to the instructions, contact the manufacturer or his responsible representative immediately.

When used as intended, these products are more or less subject to wear, depending on the intensity of use. This kind of wear is caused by technical factors and cannot be avoided. Therefore carry out a thorough visual examination and a function test prior to each application of these products. Carry out all plug connections only with the compatible articles and plug in the plug completely.

Products must be checked for irregularities. Paying attention to cracks, fractures and the occurrence of corrosion. Carry out a function check in accordance with the intended application of the product.

Essential test for bipolar electrodes is among others:

- Intact tips and blades

Defective products must not be used and must have undergone the complete preparatory treatment process again before being returned.

9 Sterilization

Prior to sterilization, products must undergo cleaning and disinfection, be rinsed off without residue using demineralized water and subsequently dried. HEBUmedical recommends using a validated steam sterilization process (e.g. sterilizer in compliance with EN 285 and validated in accordance with DIN EN ISO 17665-1).

On using the fractionated vacuum method, **sterilization** must be performed with at least **134°C (USA 132° C)** with a **minimum dwell period of 3 minutes**. Vacuum drying must then be carried out for at least 20 minutes.

 Parameters used for the validation of steam sterilization	
Prevacuum	3 times
Sterilization temperature	132 °C
Sterilization time	1,5 minutes (half cycle method)
Drying time	20 minutes

The vapour must be free of ingredients, recommended limiting values of feed water and vapour condensate are determined through EN 285.

Other sterilization processes are compatible but not validated from HEBUmedical.

When loading, observe the recommended total weight. After the sterilization, check the sterile product packaging for damage, and inspect the sterilization indicators.

The bipolar electrodes must not be sterilized in hot air sterilizers.

9.1 Packaging

Compliant packaging of products for sterilization in line with ISO 11607. Packaging used must be suitable for the instruments and protect them from microbiological contamination during storage. The seal must not be under tension. HEBUmedical recommends container or hospital common sterilization paper/film packagings for sterilized packaging.

During validation the instruments were packaged in hospital common sterilization packagings (paper/film packagings) and steam sterilized.

10 Lifetime

The steam sterilization procedure was validated by laboratory tests. The products were sterile validated at a pre-vacuum of at least 5min duration and a temperature of 134°C for a lifetime of 50 cycles.

You can continue to use the instruments at your own responsibility over this cycle value if the tests described in chapter 7 have been successfully completed.

11 Storage

Storage of the sterilized instruments in a dry, clean and dust-free environment. Products must be protected from direct sunlight in storage.

12 Warranty/Repair

Our products are made of high-quality materials and are carefully checked before delivery. However, even when used as intended, they are subject to wear to a greater or lesser extent depending on the intensity of use.

This wear is due to technical reasons and is unavoidable.

However, if faults occur that are not related to wear, please contact our customer service department. Defective products may no longer be used and must have undergone the entire preparatory process before being returned.

13 Servicing and manufacturer address

Should you require the instructions for use in paper form, please use the contact details below. The instructions for use in paper form will be made available to you within seven calendar days of receipt of the request.

Alternatively, you can print out the electronic instructions for use yourself.

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**Our vast range of products covers over 10 000 various instruments.
They can be found in Internet or request our catalog.**

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HEBUmedical GmbH